


Hans-Joachim Hübschmann

 WILEY-VCH

Handbook of GC/MS

Fundamentals and Applications

Second,
Completely
Revised and
Updated Edition

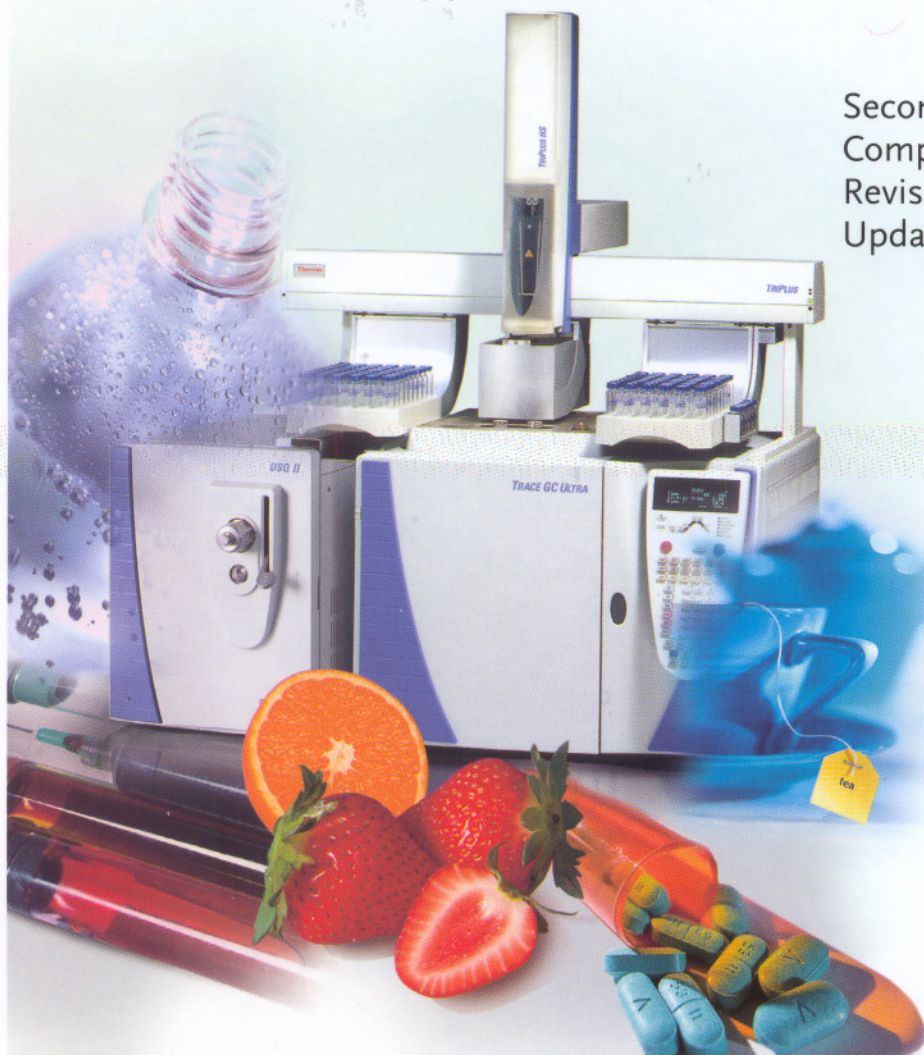

Hans-Joachim Hübschmann



Instituto Nacional
de Tecnología Industrial

Extensión y Desarrollo
División Biblioteca

13 NOV 2009

13 328 9

Handbook of GC/MS

Fundamentals and Applications

Second, Completely Revised and Updated Edition



WILEY-
VCH

WILEY-VCH Verlag GmbH & Co. KGaA

Contents

1	Introduction	1
2	Fundamentals	7
2.1	Sample Preparation	7
2.1.1	Solid Phase Extraction	10
2.1.1.1	Solid Phase Microextraction	12
2.1.2	Supercritical Fluid Extraction	15
2.1.3	Pressurized Fluid Extraction	26
2.1.4	Online Liquid Chromatography Clean-up	29
2.1.5	Headspace Techniques	30
2.1.5.1	Static Headspace Technique	31
2.1.5.2	Dynamic Headspace Technique (Purge and Trap)	39
2.1.5.3	Headspace versus Purge and Trap	49
2.1.6	Adsorptive Enrichment and Thermodesorption	54
2.1.6.1	Sample Collection	57
2.1.6.2	Calibration	59
2.1.6.3	Desorption	60
2.1.7	Pyrolysis and Thermal Extraction	63
2.1.7.1	Foil Pyrolysis	64
2.1.7.2	Curie Point Pyrolysis	66
2.1.7.3	Thermal Extraction	68
2.2	Gas Chromatography	70
2.2.1	Fast Gas Chromatography	70
2.2.1.1	Fast Chromatography	70
2.2.1.2	Ultra Fast Chromatography	74
2.2.2	Two Dimensional Gas Chromatography	75
2.2.2.1	Heart Cutting	79
2.2.2.2	Comprehensive GC×GC	79
2.2.2.3	Modulation	83
2.2.2.4	Detection	84
2.2.2.5	Data Handling	85
2.2.2.6	Moving Capillary Stream Switching	87
2.2.3	GC/MS Sample Inlet Systems	90
2.2.3.1	Carrier Gas Regulation	91

2.2.3.2	The Microseal Septum	94
2.2.3.3	Hot Sample Injection	95
2.2.3.4	Cold Injection Systems	100
2.2.3.5	Injection Volumes	108
2.2.3.6	On-column Injection	112
2.2.3.7	Cryofocusing	116
2.2.4	Capillary Columns	118
2.2.4.1	Sample Capacity	128
2.2.4.2	Internal Diameter	129
2.2.4.3	Film Thickness	130
2.2.4.4	Column Length	131
2.2.4.5	Adjusting the Carrier Gas Flow	132
2.2.4.6	Properties of Stationary Phases	134
2.2.5	Chromatography Parameters	137
2.2.5.1	The Chromatogram and its Meaning	138
2.2.5.2	Capacity Factor k'	139
2.2.5.3	Chromatographic Resolution	140
2.2.5.4	Factors Affecting the Resolution	144
2.2.5.5	Maximum Sample Capacity	146
2.2.5.6	Peak Symmetry	146
2.2.5.7	Optimisation of Flow	147
2.2.6	Classical Detectors for GC/MS Systems	151
2.2.6.1	FID	151
2.2.6.2	NPD	153
2.2.6.3	ECD	155
2.2.6.4	PID	157
2.2.6.5	ELCD	159
2.2.6.6	FPD	161
2.2.6.7	PDD	162
2.2.6.8	Connection of Classical Detectors Parallel to the Mass Spectrometer	164
2.3	Mass Spectrometry	166
2.3.1	Resolving Power and Resolution in Mass Spectrometry	167
2.3.1.1	High Resolution	174
2.3.1.2	Unit Mass Resolution	178
2.3.1.3	High and Low Resolution in the Case of Dioxin Analysis	181
2.3.2	Time-of-Flight Analyser	183
2.3.3	Isotope Ratio Monitoring GC/MS	188
2.3.4	Ionisation Procedures	206
2.3.4.1	Electron Impact Ionisation	206
2.3.4.2	Chemical Ionisation	212
2.3.5	Measuring Techniques in GC/MS	231
2.3.5.1	Detection of the Complete Spectrum (Full Scan)	231
2.3.5.2	Recording Individual Masses (SIM/MID)	233
2.3.5.3	High Resolution Accurate Mass MID Data Acquisition	246
2.3.6	MS/MS – Tandem Mass Spectrometry	250
2.3.7	Mass Calibration	261

2.4	Special Aspects of GC/MS Coupling	269
2.4.1	Vacuum Systems	269
2.4.2	GC/MS Interface Solutions	274
2.4.2.1	Open Split Coupling	274
2.4.2.2	Direct Coupling	276
2.4.2.3	Separator Techniques	277
	References for Chapter 2	278
3	Evaluation of GC/MS Analyses	293
3.1	Display of Chromatograms	293
3.1.1	Total Ion Current Chromatograms	294
3.1.2	Mass Chromatograms	295
3.2	Substance Identification	297
3.2.1	Extraction of Mass Spectra	297
3.2.2	The Retention Index	309
3.2.3	Libraries of Mass Spectra	313
3.2.3.1	Universal Mass Spectral Libraries	314
3.2.3.2	Application Libraries of Mass Spectra	317
3.2.4	Library Search Procedures	320
3.2.4.1	The INCOS/NIST Search Procedure	321
3.2.4.2	The PBM Search Procedure	328
3.2.4.3	The SISCOM Procedure	331
3.2.5	Interpretation of Mass Spectra	335
3.2.5.1	Isotope Patterns	337
3.2.5.2	Fragmentation and Rearrangement Reactions	343
3.2.5.3	DMOX Derivatives for Location of Double Bond Positions	350
3.2.6	Mass Spectroscopic Features of Selected Substance Classes	351
3.2.6.1	Volatile Halogenated Hydrocarbons	351
3.2.6.2	Benzene/Toluene/Ethylbenzene/Xylenes (BTEX, Alkylaromatics)	358
3.2.6.3	Polyaromatic Hydrocarbons (PAH)	358
3.2.6.4	Phenols	358
3.2.6.5	Pesticides	364
3.2.6.6	Polychlorinated Biphenyls (PCBs)	378
3.2.6.7	Polychlorinated Dioxins/Furans (PCDDs/PCDFs)	382
3.2.6.8	Drugs	383
3.2.6.9	Explosives	386
3.2.6.10	Chemical Warfare Agents	391
3.2.6.11	Brominated Flame Retardants (BFR)	394
3.3	Quantitation	395
3.3.1	Decision Limit	396
3.3.2	Limit of Detection	397
3.3.3	Limit of Quantitation	397
3.3.4	Sensitivity	399
3.3.5	The Calibration Function	399
3.3.6	Quantitation and Standardisation	401
3.3.6.1	External Standardization	401

3.3.6.2	Internal Standardisation	402
3.3.6.3	The Standard Addition Procedure	406
3.4	Frequently Occurring Impurities	407
	References for Chapter 3	415
4	Applications	421
4.1	Air Analysis According to EPA Method TO-14	421
4.2	BTEX Using Headspace GC/MS	429
4.3	Simultaneous Determination of Volatile Halogenated Hydrocarbons and BTEX	433
4.4	Static Headspace Analysis of Volatile Priority Pollutants	437
4.5	MAGIC 60 – Analysis of Volatile Organic Compounds	443
4.6	irm-GC/MS of Volatile Organic Compounds Using Purge and Trap Extraction	451
4.7	Vinyl Chloride in Drinking Water	454
4.8	Chloral Hydrate in Surface Water	458
4.9	Field Analysis of Soil Air	461
4.10	Residual Monomers and Polymerisation Additives	465
4.11	Geosmin and Methylisoborneol in Drinking Water	468
4.12	Substituted Phenols in Drinking Water	472
4.13	GC/MS/MS Target Compound Analysis of Pesticide Residues in Difficult Matrices	477
4.14	Multi-component Pesticide Analysis by MS/MS	489
4.15	Multi-method for the Determination of 239 Pesticides	498
4.16	Nitrophenol Herbicides in Water	505
4.17	Dinitrophenol Herbicides in Water	508
4.18	Hydroxybenzonitrile Herbicides in Drinking Water	514
4.19	Routine Analysis of 24 PAHs in Water and Soil	521
4.20	Fast GC Quantification of 16 EC Priority PAH Components	525
4.21	Analysis of Water Contaminants by On-line SPE-GC/MS	532
4.22	Determination of Polar Aromatic Amines by SPME	534
4.23	Congener Specific Isotope Analysis of Technical PCB Mixtures	540
4.24	Polychlorinated Biphenyls in Indoor Air	545
4.25	Confirmation Analysis of Dioxins and Dioxin-like PCBs	548
4.26	Fast GC Analysis for PCBs	554
4.27	Analysis of Brominated Flame Retardants PBDE	560
4.28	Trace Analysis of BFRs in Waste Water Using SPME-GC/MS/MS	568
4.29	Analysis of Military Waste	572
4.30	Detection of Drugs in Hair	582
4.31	Detection of Morphine Derivatives	584
4.32	Detection of Cannabis Consumption	589
4.33	Analysis of Steroid Hormones Using MS/MS	592
4.34	Determination of Prostaglandins Using MS/MS	596
4.35	Detection of Clenbuterol by CI	603
4.36	General Unknown Toxicological-chemical Analysis	607
4.37	Clofibric Acid in Aquatic Systems	611

4.38	Polycyclic Musks in Waste Water	616
4.39	Identification and Quantification of Trichothecene Mycotoxins	621
4.40	Highly Sensitive Screening and Quantification of Environmental Components Using Chemical Ionisation with Water	625
4.41	Characterization of Natural Waxes by Pyrolysis-GC/MS	629
4.42	Quantitative Determination of Acrylate Copolymer Layers	633
	References for Chapter 4	638
5	Glossary	639
	Subject Index	693

Subject Index

a

- abundance ratios 188 ff.
- accelerated solvent extraction (ASE) 26
- accuracy 407, 452
- acenaphthene 71, 119, 523
- acenaphthylene 119, 523, 627
- acephate 500
- acetamiprid 492
- acetate 474
- acetic acid-methoxyethylether 627
- acetone 8
 - ionisation potentials 158
 - PFE 27
 - SFE 21
- acetonitrile 8
- acetoxyscirpenol 622
- acetylation 472
- acetyldeoxynivalenol 622
- acetylene 158
- acidic hydrogen atoms 220
- acquisition 479, 489, 496
 - see also* data acquisition
- acridine 361
- acrylate copolymer layers 643
- activated charcoal 44
- active sample collection 57 ff.
- actuation 89
- adamsite 393
- ADBI polycyclic musks 616
- additives 465
- adduct formation 211, 219
- adrin 498
- adsorbents 8, 44, 54 f.
- adsorption
 - air analysis 422
 - chromatography 8
 - enrichment 54 ff.
 - headspace technique 36
 - SFE 18
 - soil air 461
 - tube calibration 61
- aerosol 102
- aggregate states 17
- agilent 120, 124
- agrochemicals 319
- AHTN 616
- air analysis 54, 421 f., 545
- alachlor 479, 501
- alcohols 216, 345, 628
- aldehyde compound classes 628
- aldrin 366, 479, 501
- Alexander Yarkov libraries 318
- aliphatic amines 215
- aliquots 7, 31
- alkanes
 - compound class 628
 - discrimination 95
 - Kovats index 309
 - proton affinities 215
 - pyrolysis 67
- alkenes 345 f.
- alkylamines 216
- alkylaromatics 345, 358
- alkylbis(trifluoromethyl) phosphine sulfides 310
- alkylphenols 472
- alkyls/derivatives 237, 313
- allethrin 492, 502
- allidochlor 499
- Alltech capillary columns 121, 125
- aluminium oxide columns 8
- amberlite 55
 - see also* XAD
- ametryn 480, 491, 501

- amines
 - cleavage 344
 - proton affinities 216
 - SPME 534
 - amino columns 8
 - aminobiphenyl 539
 - aminocarb 499, 501
 - aminotoluenes 539
 - amitriptylines 384, 610
 - ammonia 218, 227
 - amphetamines 14, 383, 582 ff.
 - anabolic steroids 3, 592
 - analysis samples *see* samples
 - analytical pyrolysis 65
 - androgens 319
 - aniline 207, 216, 539
 - anilofos 492
 - anionic samples 8
 - anthracene
 - capillary columns 119
 - *m/z* values 237
 - mass spectra 360
 - PAHs calibration 523
 - response gains 627
 - retention time 71
 - anthropogenic organic contaminants 540
 - apparatus 12, 42
 - see also* instrumental techniques
 - appearance potential 209
 - applications 421–638
 - adsorption materials 56
 - capillary columns 124
 - libraries 317 f.
 - MCSS 89
 - aqueous buffer 9
 - arachidonic acid derivatives
 - (prostaglandins) 596
 - aramite 503
 - argon
 - charge exchange 219
 - ionisation detector 5
 - ionisation potentials 158
 - aroclor 543 f., 556 f.
 - aromatic amines 534
 - aromatic hydrocarbons
 - capillary columns 119
 - chlorinated (HNU) 157 f.
 - halogenated 351
 - polycyclic 521 ff.
 - stationary phases 134
 - aromatics compound classes 628
 - aspon 502
 - asymmetric chromatogram peaks 146
 - at notation (IRMS) 191
 - atmospheric air 205
 - atrazine 312, 368, 479, 501
 - attomole region 1
 - automated mass spectra deconvolution
 - identification system (AMDIS) 305
 - automation
 - BTEX 429
 - headspace technique 53
 - mixture correction 334
 - purge and trap systems 54
 - sample inlet 90
 - sample preparation 7
 - toxicological-chemical analysis 607
 - water analysis 532
 - AutoQuan 609
 - autosampler
 - air analysis 423
 - EI/CI 231
 - volatile priority pollutants 437
 - MAGIC 60 analysis 444
 - azacozazole 492
 - azide hydrogen methylation 374
 - aziphos 480, 505
 - azulene 157, 237
- b**
- back pressure regulation 93
 - background subtraction 304
 - baking out phase 43 ff.
 - BARBs drug screening 585
 - base ion 211
 - base/neutrals and acids (BNAs) 28
 - bases 217
 - basics 7–292
 - beeswax 629
 - benalaxyl 480, 504
 - benchtop GC/MS instruments 2
 - benchtop/PBM spectral comparison 328
 - bendiocarb 370, 479, 499 f.
 - benfluralin 491, 500
 - benodanil 504
 - benoxacor 491
 - benzene
 - air analysis 428
 - breakthrough volumes 57
 - charge exchange 219

- detection limits 441
- environmental components 626
- headspace technique 36, 429
- ionisation potentials 158, 207
- irm-GC/MS 452
- ISTD 404
- m/z values 237
- MAGIC 60 analysis 447
- mass spectra 4, 355 ff.
- purge and trap/GC/MS systems 51
- recoveries 432
- benzene hexachloride (BHC)
 - acquisition parameters 479
 - multicomponent pesticide analysis 491
 - retention times 500 f.
 - water analysis 532
- benzenesulfonic acid columns 8
- benzoanthracene/fluoranthene
 - capillary columns 119
 - internal standards 527
 - PAHs calibration 523
 - response gains 627
 - retention time 71
- benzofluorine 527
- benzol 425
- benzoperylene 71, 119, 523, 527
- benzopyrene
 - capillary columns 119
 - internal standards 527
 - mass spectra 361
 - PAHs calibration 523 ff.
 - response gains 627
 - retention time 71
- benzoximate 499
- benzoylprop-ethyl 504
- benzyl cleavage 345
- benzylchlorid 425
- BFRs analysis 560
- BHG 479
- Biemann-Watson separator 277
- bifenox 504
- bifenthrin 504
- binapacryl 508 f.
- biochemical processes 3
- biomarkers libraries 317
- bioterrorism 621
- biperiden 610
- biphenyls
 - m/z values 237
 - PAHs calibration 523
 - polychlorinated 545
 - response gains 627
- bipolar fibres 14
- biscyanopropyls 136
- bitertanol 493
- blank values 53
- blood analysis 124, 553, 586
- boiling points 100
 - capillary columns 134
 - cryofocusing 117
 - injector system 112
 - PTV 103
 - retention index 313
 - solvents 107
- bond breaking 206
- bornaprin 610
- branched carbon chains 345
- breakthrough temperatures 117
- breakthrough volumes (BTv) 54, 57
- bromacil 491, 503
- bromazepam 610
- brominated aromatic compounds 394
- brominated diphenyl ethers (BDEs) 560
- brominated flame retardants 560
- bromine isotope pattern 338, 341
- bromkal 394
- bromobenzene 36, 405, 447
- bromobutane 158
- bromobutide 491
- bromochloromethane 424, 447
 - detection limits 441
 - MAGIC 60 analysis 447
 - purge and trap/GC/MS systems 51
 - recoveries 432
- bromofluorobenzene
 - air analysis 424
 - ISTD 404
 - MAGIC 60 analysis 448
- bromoform
 - breakthrough volumes 57
 - detection limits 441
 - ISTD 404
 - MAGIC 60 analysis 447
 - purge and trap/GC/MS systems 51
 - recoveries 432
- bromophenols 472
- bromophos 502
- bromopropylate 492, 504
- bromostyrene 467
- bromoxynil 506, 514

BTEX

- headspace GC/MS 429
- MAGIC 60 analysis 443
- mass spectra 358f.
- PID 158
- purge and trap systems 45
- soil air 461
- volatile halogenated hydrocarbons 433f.

Budzikiewicz charge concept 206

bufencarb 500

bulk sample isotope analysis (BSIA) 195

bupirimate 492, 504

butachlor 492, 503

butamifos 492

butane 207

butanol 35

butralin 502

butyl acetate 35, 158

butyl columns 9

butyl rubber septa 409

butylamine 158

butylate 499

butylbenzene 447

butyleneglycol-acetate 627

butylglycol 627

BZE 585

c

C¹² reference mass spectra 337

¹³C variations 193

cafenstrole 493

caffeine 610, 625

calibration

- acrylate copolymer layers 647
- air analysis 423
- BTEX 431
- chloral hydrate 459
- clenbuterol 606
- drug screening 586
- geosmin/methylisoborneol 472
- mass spectroscopy 399
- MID 247
- MS/MS 261
- multicomponent pesticide analysis 491
- PAHs analysis 523
- polychlorinated biphenyls 546
- sampling 59
- standard addition 405
- static headspace analysis 442
- toxicological-chemical analysis 609

- trichothecene mycotoxins 624

- unit 60

- vinyl chloride determination 454

calibration correction method (CCM) 173

cannabis 582, 589

capacity factor 34, 139ff., 149

capillary columns 7, 118ff.

- fast GC analysis 555

- filigree 112

- PAHs 528

- prostaglandins 598

capillary gas chromatography 1, 252

captafol 500

captan 480f., 500, 503

carbamates 370, 498

carbamazepine 610

carbenium 345

carbetamide 503

carbon 338, 345

carbon dioxide

- elution profiles 196

- ionisation potentials 207

- SFE 16ff.

carbon disulfide 441

carbon isotopic analysis 540

carbon monoxide

- charge exchange 219

- cleavage 345

- resolution 178

carbon tetrachloride

- breakthrough volumes 57

- detection limits 441

- ISTD 404

- MAGIC 60 analysis 448

- purge and trap/GC/MS systems 51

- recoveries 432

carbonyl compounds 217, 345

carbopak 45

carbophenothion 504

carborane 554, 576

carbosieve 45

carbotrap 55, 61

carbowax 14, 136

carboxin 503

carboxylic acid columns 8

carcinogenic compounds 525, 572

carnauba wax 629

carrier gas

- contaminants 408

- GC/MS coupling 270

- irm-GC/MS 451
- multi-residue method 498
- PDD 163
- regulation 91, 132
- streams 3
- velocity 130, 149
- cartridges 8
- cathode damage 274
- cationic samples 8
- C-C bonds 63, 152
- C-H bonds 152
- charge exchange 219f.
- chemical analysis 1
- chemical attack 93
- Chemical Concepts Libraries 317, 335
- chemical ionisation (CI) 6, 206, 212f.
- clenbuterol 603
- drinking water 472
- environmental components 625
- military waste 572
- PAHs analysis 521
- trichothecene mycotoxins 621
- trinitrobenzenes 579
- warfare agents 391
- water analysis 508ff.
- chemical noise 251, 303
- chemical reactions 335
- chemical warfare agents 371, 391f.
- chemometric processing 86
- chlopen 542
- chloral hydrate 458
- chlorbenside 502
- chlorbenzilate 492
- chlorbromuron 503
- chlorbufam 501
- chlordan 479, 482, 502
- chlordiazepoxide 610
- chlordimeform 500
- chlorfenon 503
- chlorfenvinphos 296, 502
- chlorflurecol-methyl 503
- chloridazon 505
- chlorinated dibenzodioxins/furans 552
- chlorinated hydrocarbons 364
- chlorination 157
- chlorine 338f.
- chlormephos 499
- chlormethan 425
- chlormezanone 610
- chloroacetophenone (CN) 393
- chloroaniline 539
- chlorobenzene
 - air analysis 424
 - breakthrough volumes 57
 - detection limits 441
 - ionisation potentials 207
 - ISTD 404
 - MAGIC 60 analysis 447
 - purge and trap/GC/MS systems 51
 - recoveries 432
 - soil air 464
- chlorobenzilate 503
- chlorobenzylidenemalonitrile (CS) 393
- chlorobromopropane 448
- chlorobutane 158
- chlorofes 542
- chlorofluoroaniline 539
- chloroform 8
 - air analysis 425
 - boiling points 107
 - breakthrough volumes 57
 - chloral hydrate 458
 - detection limits 441
 - irm-GC/MS 452
 - ISTD 404
 - MAGIC 60 analysis 447
 - mass spectra 352
 - purge and trap/GC/MS systems 51
 - recoveries 432
- chloroiodobenzenes 539
- chloroneb 500
- chlorophenol 362, 472
- chloropropane 207
- chloropropylate 503
- chlorothalonil 501
- chlorotoluene 441, 447
- chloroxynil 514
- chlorpromazine 610
- chlorpropham 500
- chlorprothixene 610
- chlorpyrifos 491, 501f.
- chlorthal-dimethyl 491, 502
- chlorthiamid 499, 502
- chlorthion 502
- chlorthiophos 504
- chlozolate 503
- cholesterol 67, 642
- cholesterylstearate 640
- chromatograms 3, 138, 293
 - PCB/Ugilec T 243

- PFB esters 614
- polychlorinated biphenyls 547
- sample trace 63
- volatile substances 47
- chrysene
 - capillary columns 119
 - internal standards 523, 527
 - response gains 627
 - retention time 71
- cigarette smoke extract 87
- clean-up 8, 405, 498
- cleavage
 - C-C bonds 63
 - ionisation procedures 208
 - library search 321
 - mass spectra 344 ff., 383
- clenbuterol 603
- clofibric acid 611
- clogging 93
- clomazone 500
- clomipramine 610
- Cl-organic compounds classes 628
- clozapine 610
- clusters 224, 321
- C-N structure 154
- coating 14, 643
- coaxial heater restrictor 20
- cocaine 385, 582, 610
- codeine 384, 584 ff., 610
- co-elution
 - chromatograms 298
 - library search 332
 - mass spectra 307
 - prostaglandins 596
- coking plant soil 322
- cola sample 624
- cold injection system 21, 100 f., 106
- cold needle 96
- cold trap 19, 62
- collision cell 255
- collision energy ramp 260
- collision-induced dissociation (CID) 253 ff.
- columns 1 ff., 7
 - bleed 118, 131, 408
 - BTEX 429, 433
 - diameters 129
 - fast gas chromatography 70, 74
 - heating technique 74
 - irm-GC/MS 451
 - length 131
 - MAGIC 60 analysis 445
 - military waste 572
 - see also* capillary columns
 - SPE 8
 - switching schematics 90
 - thick film 421
- complete spectrum detection 231
 - see also* full scan
- component perception 305
- compound classes 628
- compound derivatisation 192
- compound specific isotope analysis (CSIA) 195 ff., 451, 540
- comprehensive gas chromatography 1, 75 f.
- compression ratio 271
- concentration ranges 7
- concentrators 423, 433, 444
- concurrent solvent recondensation (CSR) 92, 100
- confirmation
 - dioxins/dioxin-like PCBs 548
 - MS/MS 489, 496
 - target compounds 477
- congener specific isotope analysis 540
- construction
 - direct coupling 276
 - ECD 155
 - ELCD 160
 - FID 152
 - FPD 162
 - oil diffusion pump 273
 - packed column SPE 11
 - PDD 163
 - PID 159
 - thermal extractor 69
 - turbo molecular pump 271
- contamination
 - air analysis 424
 - clofibric acid 611
 - congener specific isotope analysis 540
 - drinking water 468
 - frequently occurring 408 ff.
 - indoor air 545
 - irm-GC/MS 451
 - open split coupling 274
 - slow on-column injections 114
- continuous pyrolyser 65
- conversion isotope analysis 197
- conversion rates 157
- cooking oils 51

- coplanar PCBs 89
 - copolymer layers 643
 - copper 409
 - coronene 361
 - correlations, SISCOTM 334
 - cotinine 610
 - cough mixture 3
 - coumaphos 505
 - coupling
 - detectors/MS 164
 - gas chromatography/mass spectrometry 3, 270 ff.
 - purge and trap/GC/MS systems 47
 - pyrolysers/GC/MS systems 67, 70
 - cresol 362
 - critical temperatures 18
 - cross-contamination 36, 52
 - crotoxyphos 503
 - crude oil steroid biomarkers 196
 - crufomate 502
 - cryo-enrichment 108
 - cryofocusing 54, 116 f.
 - air analysis 421 ff.
 - calibration 62
 - cryogenic zone 84
 - cumenes 57
 - Curie point pyrolysis 64 ff.
 - cutoff point 164
 - cyanazine 480, 491, 503
 - cyanide/cyanate radicals 154
 - ciano columns 8
 - cyanophos 491, 501
 - cyanopropylphenyl 136
 - cycle time 249
 - cycloate 500
 - cyclodex columns 121
 - cyclohexane 34 f., 107
 - cyclohexyl columns 8
 - cyclopentapyrene 527
 - cyfluthrin 505
 - cyhalofop-buthyl 493
 - cyhalothrin 493
 - cypermethrin 493, 505
 - cyprazine 501
- d**
- data acquisition 479
 - air analysis 425
 - MDGC 85 ff.
 - MID 246 ff., 260 ff., 531
 - MS/MS 260, 489, 496
 - data file formats 308, 315
 - DDD 365, 479, 503 f.
 - DDE/DDT
 - acquisition parameters 479
 - mass spectra 366
 - minimum resolution 181
 - retention times 503
 - dead time 139 f., 147 f.
 - deca-BDE analysis 560
 - decachlorobiphenyl 381
 - decision limit 396
 - decomposition 63
 - explosives 386
 - injector system 111 ff.
 - mass spectra 370
 - deconvolution 84 ff., 305 ff., 331
 - degassing products 408
 - degradation
 - irm-GC/MS 451
 - military waste 572
 - multi-residue method 498
 - polycyclic musks 616
 - delor 542
 - delta notation/variations 189 ff.
 - delta¹³C values 540
 - deltaC/N determinations 198
 - deltamethrin 493, 505
 - demeton 479, 500
 - densities 17
 - deoxynivalenol 622
 - depletion 192
 - derivatisation 13
 - cannabis 589
 - clofibric acid 613
 - morphine 584
 - prostaglandins 596
 - SPME 534
 - water analysis 506
 - descriptor 249
 - desethylatrazin 479, 500
 - designer drugs libraries 318
 - desmethryn 480
 - desmethylclomipramine 610
 - desmethylmedazepam 610
 - desmetryn 486, 501
 - desorption
 - calibration unit 61
 - drinking water 454
 - irm-GC/MS 451

- MAGIC 60 analysis 445
- purge and trap systems 41 ff.
- detection limits (LOD) 14
- air analysis 425
- mass spectroscopy 397
- multi-residue method 498
- volatile priority pollutants 438
- water irm-GC/MS 452
- detection techniques 1
- GC/MS systems 151
- MDGC 84
- nitroaromatics 573 f.
- detoxification 607
- deuterated compounds 404
- deuterated standards 521
- DGC 540
- diacetone alcohol 410
- diacetoxyscirpenol 622
- diagnostic toxicological-chemical analysis 607
- diagram
- ion trap analyser 180
- quadrupole mass spectrometer 178
- see also* construction, schematics, principles
- dialifos 505
- di-allate 500
- diameter, capillary columns 129, 129
- diaminotoluene (DAT) 389, 539
- diazepam 610
- diazinon 188, 491, 497 ff.
- diazotation 537
- dibenzoanthracene
- capillary columns 119
- internal standards 527
- PAHs calibration 523
- response gains 627
- retention time 71
- dibenzodioxins/furans 552
- dibenzopyrene 527
- dibromochloroethane 432
- dibromochloromethane
- detection limits 441
- ISTD 404
- MAGIC 60 analysis 447
- mass spectra 352
- purge and trap/GC/MS systems 51
- dibromochlorophenol 363
- dibromochloropropane 447
- dibromoethane 354, 447
- dibutyl phthalate 413
- dichlobenil 499
- dichlofenthion 491, 501
- dichlofluanid 502
- dichlormid 499
- dichloroaniline 539
- dichlorobenzene
- detection limits 441
- internal standards 523
- ISTD 404
- MAGIC 60 analysis 447 f.
- mass spectra 355
- dichlorobiphenyl 379
- dichlorobromomethane 405
- dichlorodifluoromethane 448
- dichloroethane
- breakthrough volumes 57
- detection limits 441
- MAGIC 60 analysis 448
- purge and trap/GC/MS systems 51
- recoveries 432
- dichloroethylene
- calibration file 446
- irm-GC/MS 452
- MAGIC 60 analysis 448
- recoveries 432
- static headspace analysis 440
- dichloriodobenzene 539
- dichloromethane 8, 16
- dichloromethane
- boiling points 107
- detection limits 441
- ISTD 404
- MAGIC 60 analysis 448
- mass spectra 352
- PAHs calibration 523
- PFE 27
- purge and trap/GC/MS systems 51
- recoveries 432
- dichlorophenol 363
- dichlorophenyl acetate 363
- dichloropropane 57, 432, 448
- dichloropropene 354, 441
- dichloropropylene 448, 613
- dichlorvos 479, 499
- diclofop-methyl 492, 504
- dicloran 491, 501
- dicofol, 502
- dicrotophos 500
- dieldrin 367, 479, 503

- diethofencarb 491
 - diethyleneglycols 627
 - diethylether 107
 - diffusion 139f., 147f.
 - sampling tubes 57
 - SFE 17
 - difluorobenzene 424
 - digital image processing 87
 - dihalogenated hydroxybenzonitrile 514
 - dihydrocodeine/morphine-TMS 587
 - diiodo(nitro)toluene 539
 - dilution 7, 112
 - dimepiperate 492
 - dimetachlor 501
 - dimethenamid 491
 - dimethoate 479, 501
 - dimethyl 121
 - dimethyl disulfide (DMDS) 350
 - dimethylformamide 158
 - dimethylnaphthalene 359, 523, 627
 - dimethylnitrophenol 506
 - dimethyloxazoline (DMOX) derivatives 320, 350
 - dimethyl-polysiloxane 134
 - dimpylate 479
 - dinitramine 501
 - dinitrobenzene (DNB) 390ff.
 - dinitrophenol herbicides 508f.
 - dinitrotoluene (DNT) 388
 - dinobuton 508f.
 - dinoseb/terb herbicides 508ff.
 - dioctyl phthalate 413
 - diol columns 8
 - dioxan 34
 - dioxathion 501
 - dioxins 3
 - capillary columns 124
 - confirmation analysis 548
 - ion trap analyser 181
 - PCBs 548, 554
 - PFE 28
 - polychlorinated 382
 - diphenamid 491, 502
 - diphenylamine 500
 - diphenyl-dimethyl-polysiloxane 134
 - diphenylhydramine 610
 - direct GC/MS coupling 276ff.
 - discharge gas 163
 - discrimination 95
 - dispersion 145
 - dissociation 155, 191
 - distribution constant 14
 - disubstituted alkanes 215
 - disubstituted aromatics 313
 - disulfoton 188, 306, 501
 - dithiophosphoric acid esters 371
 - diuron 375
 - divinylbenzene/carboxen/polydimethylsiloxane (DVB/CAR/PDMS) 14
 - DNOC structures 509
 - doping analysis, steroid hormones 592
 - double bonds 346, 350
 - double focusing magnetic sector instrument 1
 - doxepin 610
 - doxylamine 610
 - drinking water, *see also* water
 - analysis 405
 - geosmin/methylisoborneol 468
 - herbicides 508
 - MAGIC 60 analysis 443
 - SPME 534
 - substituted phenols 472
 - vinyl chloride 454
 - droplet formation 102
 - drugs
 - acrylate copolymer layers 643
 - hair analysis 582
 - libraries 319
 - mass spectra 383
 - morphine derivatives 584
 - toxicological-chemical analysis 607
 - drying 15, 42ff.
 - dual cryo jet modulator 84
 - duty cycles 232, 250
 - dwell times 232f., 247
 - dynamic headspace technique (purge trap) 39
 - dynamic regions 521
- e**
- Eddy diffusion 139, 146
 - edifenphos 492, 504
 - effective plates 150
 - electrolytical conductivity detector (ELCD) 159f.
 - electron beam ionisation 335
 - electron capture 155, 222
 - electron capture detector (ECD) 12, 51–157
 - military waste 572f.

- target compound analysis 477
 - electron impact ionisation (EI) 4, 206
 - environmental components 625
 - MS/MS 252
 - perfluorophenanthrene 268
 - trichothecene mycotoxins 621
 - trinitrobenzenes 579
 - electron loss 344
 - electron migration 206
 - electronic flow control units (EPC) 91
 - electrospray ionisation (ESI) 183
 - elemental analysers (EA) 194
 - elution/eluates 8
 - air analysis 424
 - capillary columns 133
 - chromatograms 139, 293
 - internal standard 403
 - IRMS 196
 - MAGIC 60 analysis 443
 - SPE 8
 - emergency analysis 607
 - endosulfan 479, 492, 502 ff.
 - endrin 367, 479, 503
 - enrichment 11
 - environmental components 625
 - environmental congener specific isotope analysis 540
 - Environmental Protection Agency (EPA)
 - method 16, 46
 - air analysis 421 ff.
 - MAGIC 60 analysis 443
 - PAHs analysis 521
 - enzymatic reactions 192
 - ephedrine 383
 - EPN 504
 - EPTC 499
 - equilibrium fractionation 191
 - erbon 503
 - errors, standard addition 407
 - ester compound class 628
 - etacnazole 504
 - ethalfuralin 491, 500
 - ethane 158
 - ethers 344, 628
 - ethion 504
 - ethofumesate 502
 - ethoprophos 479, 491, 500
 - ethylacetate 8
 - ethylacetate
 - boiling points 107
 - partition coefficients 34
 - SFE 21
 - ethylacrylate layers 643
 - ethylan 503
 - ethylbenzene
 - breakthrough volumes 57
 - detection limits 441
 - headspace technique 36
 - ISTD 404
 - MAGIC 60 analysis 448
 - mass spectra 356 ff.
 - purge and trap/GC/MS systems 51
 - recoveries 432
 - ethylchloride 447
 - ethylene 158, 176
 - ethyleneglycols 627
 - ethyl-methylbenzene 357
 - ethylnaphthalene 359
 - etoloxamine 610
 - etoxazole 492
 - etridiazole 499
 - etrimfos 480, 501
 - EU regulation 525, 534
 - evaluation, GC/MS analyses 293–420
 - evaporation rate 100
 - exact masses 338
 - excitation voltage 479 f.
 - excited radicals 151
 - expansion volumes 110
 - explosives 386
 - exponential mass calibration 262
 - external standardisation 401, 454, 459
 - extraction 7
 - mass spectra 297
 - nitroaromatics 572
 - PFE 27
 - SFE 16
 - standard addition 405
- f**
- Faraday cup arrangement 203
 - fast gas chromatography 70 ff., 525, 554
 - fast on-column injections 114
 - fate 540
 - fatty acid methyl ester (FAME) 124, 201
 - fatty acids 320, 640
 - FC43 spectrum 264 ff.
 - fenamiphos 503
 - fenarimol 480, 505
 - fenchlorphos 501

- fenitrothion 502
- fenoxanil 492
- fenpropathrin 492, 504
- fenson 502
- fensulfothion 504
- fenthion 502
- fenvalerate 493, 505
- ferromagnetic alloys 66
- fibres types 14
- field analysis 461
- filigree capillary columns 112
- filled needle hot sample injection 96
- film thickness 118, 130f.
- fipronil 491
- fish samples 481, 620
- FIT values 324f.
- flame ionisation detector (FID) 151, 294
- flame retardants 560
- flamephotometric detector (FPD) 161
- flameproofing agents 394
- flamprop-isopropyl 504
- flamprop-methyl 492, 503
- flash pyrolysis 65
- flavour components 202, 525
- flight times 211
- Florasil columns 8
- flow diagram 19
- flow divider 165
- flow rates 91ff.
 - capillary columns 130
 - chromatography parameters 138
 - fast gas chromatography 72
 - optimisation 147
 - PDD 164
- flow switching 84, 87
- fluacrypyrim 492
- fluazifob-buthyl 480
- fluchloralin 501
- flumetralin 503
- flumioxazin 493
- fluoranthene
 - capillary columns 119
 - PAHs calibration 523
 - response gains 627
 - retention time 71
- fluorene
 - capillary columns 119
 - isotope frequencies 338
 - m/z values 237
 - PAHs calibration 523
 - response gains 627
 - retention time 71
- fluorobenzene 405, 448
- fluorodifen 504
- fluorochloridone 503
- flushing air analysis 424
- flutlanil 492
- fly ash 22f.
- foil pyrolysis 62ff.
- folpet 499, 503
- fonofos 500
- food
 - chain 3
 - MS/MS 489
 - PAHs quantification 525
 - trichothecene mycotoxins 621
 - volatiles libraries 318
- formaldehyde 14
- fortifying 582
- forward pressure regulation 91
- fosthiazate 491
- fractionated pyrolysis 65
- fractionation 191
 - see also* main fragment m/z values
- fragmentation
 - electron impact ionisation 209
 - GC/MS systems 314
 - mass spectra 335ff., 343f., 383
 - MS 211
 - phosphoric acid esters 371
 - prostaglandins 599
 - protonation 218
- freezing/freeze drying 15ff.
- freons 16, 425
- frit spargers 40
- fruits 498
- FT Orbitrap mass analysers 174
- fthalide 491
- fuels 89
- full scan 231f., 468
- full width at half maximum (FWHM) 167
- fundamentals 7–292
- furans 3
 - capillary columns 124
 - confirmation analysis 553
 - PFE 28
 - polychlorinated 382
- fusarenon 622
- fused silica capillary columns 2–8, 70, 90, 120

g

- gas chromatography (GC) 1, 70 ff.
- air analysis 423
- see also* GC/MS
- gas extraction techniques 32
- gas lines 424
- gas liquid chromatography (GLC) 4
- gas phase 17
- gas phase acidities 221
- Gaussian curves 138
- GC/MS techniques 1 ff., 270 ff.
- air analysis 423
- BTEX determination 433
- prostaglandins 596
- steroid hormones 592
- trichothecene mycotoxins 621
- GC/MS/MS techniques 254
- PBBs 568
- pesticides chromatograms 495
- target compound analysis 477
- geochemicals mass spectra libraries 317
- geosmin/methylisoborneol 468
- glass columns 8
- glass dome operation 87
- glass PAT apparatus 42
- glycols/derivatives 628
- graphitised carbon black 56
- Graphpac-D 45
- ground water 538, 611
- see also* water, waste-, surface-, drinking-

h

- ²H variations 194
- hair analysis 582
- Halfenprox 493
- Hall detector 159
- halogen content 12
- halogenated hydrocarbons 50
- air analysis 421 ff.
- BTEX 432
- chromatograms 293, 298, 304
- internal standard 405
- MAGIC 60 analysis 443
- mass spectroscopy 351
- PID 158
- soil air 461
- haloperidol 610
- HCH 502
- HCN loss 346
- headspace technique 3, 7, 30 ff.

- BTEX 429
- PBBs 568
- vs. purge trap 49, 54
- heart cutting 75 f.
- heating rates 133
- height equivalent to theoretical plate (HETP) 141 f., 147 ff.
- helium
- GC/MS coupling 270
- ionisation potentials 158, 207
- PDD 163
- van Deemter curves 133
- viscosity 166
- water removal 199
- heptabromobiphenylene 223
- heptachlor 479, 501 f.
- heptachlorobiphenyl 181, 380
- heptane/heptene 57
- heptenophos 479, 500
- herbicides 28, 368, 505 ff.
- heroin 385, 583 ff.
- heterocycles 346
- hexabromobiphenyl (HBB) 223, 394
- hexachlorobenzene (HCB) 365, 499 f., 558, 628
- hexachlorobiphenyl 245, 380
- hexachlorobutadiene 405, 441,, 448
- hexachloroethane 441
- hexaconazole 492
- hexane 8
- boiling points 107
- ionisation potentials 158
- partition coefficients 34
- PFE 27
- hexazinon 369, 480,, 504
- hexogen (RDX) 386 f.
- HHCB 616
- high boiling compounds 2
- high frequency pyrolysis 64, 66 f.
- high resolution gas chromatography (HRGC) 90
- high resolution magnetic sector analyser 185
- high resolution mass spectrometers 168
- high resolution MID data acquisition 246 f.
- high resolution PAHs quantification 525
- high resolution quadrupole analysers 173 ff.
- high speed column conditions 84
- high temperature phase capillary columns 119
- highly sensitive environmental screening 625
- historical GC/MS technique development 4
- hit quota (library search) 330

- homeland security 621
 - homologous M-series 310
 - hormones 592, 596
 - hot sample injection 95 f.
 - hot split injector 603
 - HPLC 2, 12, 18, 584
 - HRMS 525, 548
 - H-SRM 173, 489, 496
 - see also* selected reaction monitoring (SRM)
 - HT-2 toxin 621
 - human body biochemical processes 3
 - humps 75
 - hydrazines 217
 - hydride abstraction 218
 - hydrocarbons, *see also* halogenated
 - hydrocarbons
 - aromatic/chlorinated (HNU) 158
 - background contaminants 411
 - BTEX 432
 - chlorinated 364
 - chromatograms 293, 298, 304
 - ECD 156
 - MAGIC 60 analysis 443
 - mass spectroscopy 351
 - PID 158
 - polyaromatic 2, 17, 23 f., 153
 - polynuclear aromatic 70, 119
 - pumps 273
 - soil air 461
 - stationary phases 134
 - total petroleum (TPH) 28
 - wax pyrograms 640
 - hydrocodone 587, 610
 - hydrogen
 - cleavage 346
 - FID hydrogenator 152
 - GC/MS coupling 272
 - ionisation potentials 158
 - isotope frequencies 338
 - retention index 313
 - van Deemter curves 133
 - viscosity 166
 - hydrolysis
 - cannabis detection 589
 - chloral hydrate 458
 - dinitrophenol herbicides 510
 - morphine derivatives 584
 - hydromorphone 587
 - hydroxybenzoxazole herbicides 514
 - hydroxytetrachlorodibenzofuran 181
 - hyperbolic quadrupoles 171 ff.
 - hyphenated techniques 7
- i*
- identical treatment of reference and sample
 - material (IT) 204
 - identification, trichothecene mycotoxins 621
 - image processing 87
 - imazalil 503
 - imazamethabenzmethyl 492
 - imibenconazole 493
 - iminoethers 217
 - imipramine 610
 - immunological methods 584
 - impurities 407 ff.
 - in situ derivatisation 534
 - INCOS library search 320 ff., 331 ff.
 - indenopyrene
 - capillary columns 119
 - internal standards 527
 - PAHs calibration 523
 - response gains 627
 - retention time 71
 - indoor air 545
 - see also* air, soil air
 - information content 3, 253
 - injection 89 ff.
 - air analysis 422
 - BTEX 429, 433
 - chloral hydrate 458
 - drinking water 454
 - headspace technique 37, 53
 - hot/cold 111
 - irm-GC/MS 451
 - IRMS 204
 - MAGIC 60 analysis 445
 - MS/MS 490
 - multi-residue method 498
 - PAHs analysis 521
 - pressure regulation 93 f.
 - purge and trap/GC/MS systems 47
 - pyrolysis 66
 - target compound analysis 477
 - volatile priority pollutants 437
 - volumes 108 f.
 - instrumental techniques 1
 - purge and trap systems 42
 - SPE 12
 - integrated control systems (INCOS) library 321

intensities 28, 334
 interfaces 81, 274 ff.
 interfering signals 407
 internal capillary column diameter 129
 internal standard (ISTD)
 – BTEX determination 433
 – cannabis detection 589
 – clofibric acid 611
 – MAGIC 60 analysis 443, 448
 – mass spectroscopy 402, 404
 – multi-residue method 498
 – PAHs analysis 523, 527
 – static headspace analysis 437
 interpretation, mass spectra 335 ff.
 involatile matrix 112
 iodination 534
 iodine 338
 iodobenzene/biphenyl 539
 iodobutane 158
 iodonaphthalene 539
 iodotoluenes 539
 ion abundance dependence 187
 ion beam instruments 231
 ion chromatogram *see* chromatograms
 ion cyclotron resonance 5, 174
 ion exchange chromatography (IEC) 9
 ion source stability 53
 ion storage mass spectrometers 232
 ion trace chromatogram *see* chromatograms
 ion transmission 173
 ion trap instruments 6
 – EI/CI 229–235
 – environmental components 625
 – MS/MS 252, 256
 – prostaglandins mass spectra 601
 – storage times 211
 ion types 211
 ionic samples 8
 ionisation 206
 – FID 151 f.
 – irm-GC/MS 452
 – mass spectra 335
 – MS/MS 252
 – PAHs quantification 528
 ionisation potentials 4, 158, 207 f.
 ionol (BHT) 413
 ioxynil 506 f., 514
 iprodione 504
 irm-GC/MS 451 f.
 isazofos 501

isobaric interferences 203
 isobutane 218, 227
 isobutyraldehyde 158
 isofenphos 502
 isomer separation 529
 isomeric compounds 2, 324
 isopentane 107
 isopropalin 502
 isopropanol 35
 isopropylbenzene 448
 isopropyltoluene 405, 448
 isotope dilution 403 f.
 isotope frequencies 338
 isotope patterns
 – library search 321
 – mass spectra 336 ff.
 – PCB 245
 – TCDD 182
 isotope ratio mass spectrometry (IRMS)
 88, 188 ff., 202
 isotope reference materials 205
 isotopes 4
 isotopic fractionation 191

j

jet-pulsed modulators 84

k

kanechlor 542
 Kel-F 409
 ketones 344, 628
 key fragments *see* main fragments
 kieselgur columns 8
 kinetic energy 191
 kinetic fractionation 192
 Kovats index 150

l

landfill gas 433
 lanolin 629
 large volume injection 107 f., 477
 lemon extract analysis 295
 lenacil 492
 lengths, capillary columns 131
 leptophos 504
 levallorphan 587
 levomepromazine 610
 levorphanol 587
 library search
 – air analysis 427

- mass spectra 306, 313 ff., 320, 335
 - soil air 461
 - volatile priority pollutants 437, 441
 - lidocaine 610
 - limits of detection (LOD) 1
 - air analysis 425
 - multi-residue method 498
 - volatile priority pollutants 438
 - limits of quantitation (LOQs),
see also quantitation
 - MAGIC 60 analysis 447
 - PBDE analysis 561
 - lindane
 - acquisition parameters 479
 - mass spectra 364
 - retention times 501
 - water analysis 532
 - line materials contaminants 409
 - line spectrum 210
 - linear calibration function 399
 - linear range 14
 - linear temperature programs 309
 - linuron 375, 502
 - lipids libraries 320
 - liquid extraction procedures 572
 - liquid nitrogen trap 200
 - liquid phase 17
 - liquid solid chromatography (LSC) 9
 - liquid-liquid microextraction (LLME) 437
 - local normalisation (library search) 321, 328
 - localised charge concept 206
 - lock-and-cali mass technique 248 f., 549
 - lock-mass technique 247 f.
 - lodofenphos 503
 - loop injection 422
 - low resolution quadrupole analysers 173 ff.
 - low voltage ionisation 206
 - lower application limits 51
 - LSD 585
 - Ludwig light phase rotation perforator 507
- m**
- m/z* values see main fragments
 - Macherey-Nagel capillary columns 121, 125
 - MAGIC 60 analysis 236, 443
 - magnetic sector instruments 1, 174
 - calibration 262
 - chromatograms 294
 - EI/CI 229 f.
 - flight times 211
 - TCDD isotope pattern 184
 - Magnum CI (advanced chemical ionisation) 213
 - main fragment *m/z* values 297 f., 347 ff.
see also fragments
 - contaminants 408
 - pesticides 238
 - polycondensed aromatics 237
 - MAK value list 573
 - malaoxon 502
 - malathion 372, 502
 - mass analysers 174
 - mass calibration 261
 - mass chromatograms see chromatograms
 - mass correlations 347
 - mass detection scheme 247
 - mass equilibrium 33
 - mass filter 5
 - mass instability 6
 - mass selectivity 12
 - mass spectra
 - ADBI 618
 - AHTN 618
 - carboxy-THC-TMS derivatives 590
 - chromatograms 293
 - clenbuterol 603
 - HHCB 617
 - interpretation 335 ff.
 - methandienone 593
 - PBDEs 561
 - prostaglandins 598
 - trinitrobenzenes 579
 - mass spectrometry (MS) 1 ff., 166 ff., 231
 - BTEX 433
 - isotope ratio analysis 200
 - MassLib program package 335
 - matrix 7
 - headspace technique 35
 - standard addition 405
 - target compound analysis 477
 - matrix assisted laser desorption (MALDI) 183
 - matrix exclusion 89
 - matrix noise 371
 - matrix swelling 21
 - maximum sample capacity 146
 - McLafferty rearrangement 206 f., 346 ff.
 - McLafferty statistics, library search 328
 - MCPA/B methylester 377
 - measuring techniques 231 ff.
 - meat 525

- mecoprop 514, 613
- medazepam 610
- mefenacet 493
- mefenoxam 491
- mefenpyr-diethyl 492
- Merlin microsealsepta 94
- metabolites
 - cannabis detection 589
 - mass spectra libraries 319
 - military waste 572f.
 - morphine 584f.
- metal chelates 9
- metalaxyl 501
- metamizol-desalkyl 610
- metastable ions 211
- metazachlor 502
- methadone 582, 610
- methamidophos 499
- methandienone 592
- methane
 - FID 151
 - ionisation potentials 158
 - NCI 224
 - oxidation reactor 198
 - protonation 218
- methanol 8
 - boiling points 107
 - environmental components 625
 - MAGIC 60 analysis 443
 - NCI 225
 - PFE 27
 - protonation 218
 - SFE 21
- methaqualone 610
- methidathion 480, 503
- method detection limits (MDL) 173
- methoprotrene 503
- methoximes 596
- methoxychlor 479, 504
- methyl bromide/chloride 447
- methyl esters
 - mass spectra 320, 376f.
 - EI analysis 515
 - nitrophenol herbicides 508
- methyl isothiocyanate 158
- methyl linolenate 207
- methyl trithion 503
- methylacetate 107
- methylation 505 ff., 514 ff.
- methylbutylether 452
- methylchrysene 527
- methylidiuron 375
- methylisoborneol (MIB)/geosmin 468
- methylisopropylbenzene 358
- methyllinuron 376
- methylmethacrylate copolymer layers 643
- methylmonuron 374
- methylnaphthalene 523, 627
- methylphenanthrene 523, 627
- methylphenobarbital 610
- methylsilicone phases 237
- metobromuron 501
- metolachlor 502
- metribuzin 501
- mevinphos 479, 491, 500
- mexacarbate 500f.
- microcontaminants 532
- microextraction 7
- microseal septum 94f.
- microwave assisted extraction (MAE) 26
- midazolam 610
- migration 8, 346
- military waste analysis 572
- minimum detectable quantity (MDQ) 164
- minimum resolution 176, 181
- mirex 367, 499, 504
- mixture correction 334
- mobile soil air analysis 461
- moclobemid 610
- modifiers 21ff.
- moisture control system (MCS) 43, 47, 54
- molecular ions 211, 336
- molecular sieves 56, 61
- molecular sulfur 411
- monoacetylmorphine (MAM) 583 ff.
- monoacetylmorphine-TMS 587
- monochlorobiphenyl 379
- monocrotophos 491, 501
- monolinuron 501
- mononitrobenzene (MNB) 391
- mononitrotoluene (MNT) 390
- monuron 374
- morphine 384 ff., 587
- moving capillary stream switching (MCSS) 87f., 540
- MS/MS (tandem mass spectrometry) 250 ff.
 - hair analysis 582
 - pesticides 489
 - product ion spectrum THC-TMS 591
 - prostaglandins 596

- quadrupole analysers 173
- steroid hormones 592
- trichothecene mycotoxins 621
- MTBE 107, 441
- multibed adsorption/desorption tube
 - carbotrap 300 56
- multicomponent analyses 236, 489
- multidimensional gas chromatography (MDGC) 77
- multidimensional mass spectrometry 251
- multiple headspace extraction procedure (MHE) 32
- multiple ion detection (MID) 231, 234
 - brominated flame retardants 560
 - dioxins 548
 - PAHs 531
 - PCB/Ugilec 242
 - PFB esters 614
- multiple split technique 62
- multipoint calibrations 400
- multiresidue method 489, 498
- musk perfumes 616
- myclobutanil 504
- mycotoxins 621
- n**
- ¹⁵N variations 194
- nafion 61, 200, 421
- nalorphine 587
- naloxone 610
- naphthalene 55
 - capillary columns 119
 - detection limits 441
 - internal standards 523
 - ISTD 404
 - *m/z* values 237
 - MAGIC 60 analysis 448
 - mass spectra 359
 - response gains 627
 - retention time 71
- naphthoquinone 157
- naphthyl compound class 628
- naphthylamine 539
- naphthylene 71
- natural compounds, ¹³C/²H/¹⁸O variations 193f.
- natural isotope frequencies 338
- natural waxes 629
- needle sparger 40
- needles, hot/cold/filled 95
- negative chemical ionisation (NCI) 217, 220ff.
 - mode FC43 spectrum 269
 - MS/MS 252, 258
- neoprene contaminants 409
- neosolaniol 622
- Nernst equation 138
- neutral particles loss 345
- nicotine 610
- NIST library
 - EPA/NIH 314, 321
 - PAHs analysis 523
 - volatile priority pollutants 437, 441
- nitralin 505
- nitrapyrin 500
- nitrazepam 610
- nitric oxide 219
- nitroaniline 539
- nitroaromatic compounds (NAC) 534, 573, 628
- nitrobenzene 157, 207, 539
- nitrocompounds 154, 216
- nitrofen 503
- nitrogen
 - charge exchange 219
 - compound class 628
 - ionisation potentials 158, 207
 - isotope frequencies 338
 - isotope reference materials 205
 - mass spectra 336
 - resolution 176
 - van Deemter curves 133
 - viscosity 166
- nitrogen/phosphorus detector (NPD) 154f., 477
- nitrogen/phosphorus pesticides 479
- nitrogroups 386
- nitromusk compounds 616
- nitropenta (PETN) 386ff.
- nitrophenol herbicides 505
- nitrothal-isopropyl 491, 502
- nitrotoluene 386f.
- nivalenol 621
- noise analysis 251, 305
- noise filter 321, 327
- nomenclature, pyrolysis 65
- nonachlor 479
- nonachlorobiphenyl 181, 381
- nonpolar fibres 14
- norcodeine 587
- nordazepam 610
- norflurazon 504

normal phase chromatography (NPC) 8, 9
 nortilidine 610
 nortrimipramine 610
 nortriptyline 610
 notations, IRMS 189
 nuarimol 504
 nuclide masses 4
 number of theoretical plates 149

o

¹⁸O variations 193
 octachlorobiphenyl 381
 octachlorostyrene 365
 octadecatrienoate 351
 octadecyl columns 8
 octafluoronaphthalene 222
 octogen 386
 octyl columns 8
 odour contamination 468, 472
 off-line techniques 7
 oil diffusion pump 273
 oils 51, 89
 oleic acid amide 411
 omegatron 5
 omethoate 500
 on-column injection 112
 one-point calibrations 400
 one-step jet separator 277
 on-line coupling 7, 24
 on-line high temperature conversion 200
 on-line sample combustion 197
 on-line SPE-GC/MS 532
 open split coupling, GC/MS 274
 open split interface 196 f., 423
 open tubular columns 5
 operation modes 40, 255
 opiates 582, 585
 optical resolution 179
 Orbitrap mass analysers 174
 organic acids compound class 628
 organic chemical analysis 1
 organic compounds libraries 318
 organic samples/solvents 8
 organochlorine pesticides 28, 479
 organophosphorus pesticides 28
 OV-1, purge and trap systems 44
 oven program optimisation 72
 oxadiazon 503
 oxadixyl 480, 492, 504
 oxidation reactor 198

oxidative pyrolysis 65
 oxides 215
 oxycarboxin 499, 504
 oxychlordane 502
 oxycodone 587
 oxyfluorfen 492, 504
 oxygen
 – contaminants 409
 – detector (O-FID) 152
 – ionisation potentials 207
 – isotope frequencies 338
 – proton affinities 216
 ozone problem 3

p

packed column technique 1, 71, 277 f.
 paclobutrazol 492
 Palisade Complete Mass Spectra Library
 317, 328 ff.
 panathion-methyl 502
 paracetamol 610
 paramagnetic properties 66
 paraoxon 502
 paraoxon-ethyl/methyl 372
 parathion 491, 502
 parathion-ethyl 312, 373
 parathion-methyl 373, 491
 Parlar 69 230
 partition coefficient 12
 – chromatograms 138 ff., 149
 – headspace technique 30, 33 f.
 – halogenated hydrocarbons 50
 passivation 422
 passive sample collection 57 ff.
 pattern recognition 327, 629
 patulin 621
 PBDE analysis 560, 568
 PBM (probability based match) library
 search 320, 328 f.
 PCP 585
 peak capacities 75
 peak symmetry 146
 peak width 140 ff.
 Pebulate 499
 Pee Dee Belemnite isotope reference
 materials 205
 penconazole 480, 503
 pendimethalin 502
 pentachlorobiphenyl 380
 pentachlorobiphenylene 181

- pentafluorobenzyl esters 611
- pentane 107, 207
- pentazocine-TMS 587
- pentobarbital 610
- perazine 610
- perchloro compounds classes 628
- perfluorinated compounds 261
- perfluorokerosene (PFK) 267 ff.
- perfluorophenanthrene 268
- perfluorotributylamine (FC43) 313
- perforator 505
- PerkinElmer capillary columns 121, 125
- permethrin 505
- perylene 523, 627
- pesticides 3
 - capillary columns 124
 - compound class 628
 - libraries 309, 319
 - M series chromatograms 311
 - main fragments 238
 - mass spectra 368 ff., 371 ff.
 - MCSS 89
 - multicomponent MS/MS 489, 498
 - PFE 28
 - target compound analysis 477
- Pethidine 610
- petrochemicals mass spectra libraries 317
- PFB-TMS derivatives 598
- phalates compound class 628
- pharmaceuticals 3, 319, 611
- phases
 - capillary columns 120
 - chromatograms 149
 - IRMS 191
 - purge and trap/GC/MS systems 48
 - SFE 16f.
- phenacetin 610
- phenanthrene 71, 119, 523, 627
- phenazone 610
- phenobarbital 610
- phenols
 - compound class 628
 - drinking water 472
 - ionisation potentials 158
 - mass spectra 358, 362
 - PFE 26
- phenomenex 121, 125
- phenothrin 492f.
- phenoxyalkylcarboxylic acids 376, 611
- phenthoate 492, 502
- phenylacetates 474
- phenylphenol 500
- phenyls 8, 527, 573
- phenylureas 14, 374
- phenytoin 610
- phorate 500
- phosalone 373, 504
- phosmet 504
- phosphamidon 491, 502
- phosphoric acid esters 297, 371, 391
- phosphorus 338
- phosphorus compounds 162
- phosphorus-containing substances 154
- photoionisation detector (PID) 157f., 572
- photon elimination 255
- physical properties, SFE 17
- phytane 196
- picolinyl esters 320
- piperphos 492
- pirimicarb 370, 480, 501
- pirimiphos 480, 49, 501f.
- piston pump 18
- plant protection agents 3, 295, 364
- plastic sheeting 51
- plasticisers 8, 408
- plungers 13
- poisons mass spectra libraries 319
- polar aromatic amines 534
- polar fibres 14
- polar groups 313
- polar phenols 26
- polar samples 8
- polar solvents 21
- polar stationary phases 118
- polarity 14, 468
- pollutants mass spectra libraries 319
- polyacrylate (PA) 14
- polyaromatic hydrocarbons 2, 17, 23f., 153
- polychlorinated biphenyls (PCBs)
 - capillary columns 124
 - compound class 628
 - congener specific isotope analysis 540
 - ELCD 161
 - external standardisation 401
 - indoor air 545
 - mass spectra 378
 - MCSS 89
 - PFE 28
 - processing 12
 - PTV 106

- SIM analysis 241 ff.
- polychlorinated bornane 230
- polychlorinated dioxins (PCDDs) 3, 382
 - confirmation analysis 548
- PFE 28
- SFE 22 f.
- polychlorinated furans (PCDFs) 382, 548
- polycondensed aromatics 237
- polycyclic aromatic hydrocarbons (PAHs)
 - analysis 521 ff.
 - anthracene compounds 628
 - capillary columns 119
 - fast gas chromatography 70
 - military waste 572
 - response gains 627
- polycyclic musks 616
- polydimethylsiloxane (PDMS) coating 14 f.
- polyethylacrylate-methylmethacrylate (PEAMMA) 643
- polyethylene 409
- polymerisation additives 465
- polynuclear aromatics 358
- polysiloxane capillary columns 121
- POPs 532
- Porapak 56, 61
- porous carbons 13
- porous layer open tubular (PLOT) column 454 ff.
- positive chemical ionisation (PCI) 218 ff.
 - clenbuterol 605
 - environmental components 625
 - FC43 spectrum 268
 - MS/MS 252, 270
- powdered activated carbon (PAC) 14
- precursor ions 479 f.
- precursor-product ion scan 254
- preparation (standard solutions) 404
- pressure balanced injection 37, 429
- pressure equivalence principle 15
- pressure equivalent extraction 18
- pressurized fluid extraction (PFE) 26 ff.
- pressurized solvent extraction (ASE) 526
- pretilachlor 492
- primidone 610
- principles
 - chemical ionisation 212
 - concurrent solvent recondensation 100
 - ECD 155
 - ELCD 160
 - FPD 162
 - headspace technique 31
 - high resolution adjustment 169
 - high temperature conversion 200
 - identical treatment of reference and sample material (IT) 204
 - ion trap analyser 180
 - isotope ratio monitoring 189 f.
 - magnetic sector mass spectrometer 175
 - mass spectrometry 343
 - multidimensional gas chromatography 82
 - PDD 163
 - PID 158
 - quadrupole mass spectrometer 178
 - SRM target compound analysis 259
 - TOF analyser 186
 - two dimensional gas chromatography 80
- pristane 196
- probability based match (PBM) library search 328 f.
- prochloraz 505
- procymidone 503
- profenofos 492, 503
- profluralin 501
- programmed temperature vaporiser (PTV) 100 ff., 115, 560
- prohydrojasmon 491
- promazine 610
- promecarb 370, 499 d
- promethazine 610
- prometon 500
- prometryn 480, 501
- propachlor 500
- propane 207
- propanil 491, 502
- propargite, 504
- propazine 491, 494, 501
- propene 158
- propetamphos 504
- propham 500
- propiconazole 504
- propoxur 479
- propranolol 610
- propyl phthalate 412
- propylamine 207
- propyleneglycols 627
- propyphenazones 610
- propyzamide 501
- prostaglandins 596
- prothiofos 503
- proton abstraction 220

proton affinities 215
 protonation 214, 218f.
 pulse pyrolyser 65
 pulsed discharge detector (PDD) 162f.
 pulsed positive ion negative ion chemical
 ionisation (PPINICI) 218
 purge trap techniques 3, 7
 – BTEX determination 433
 – chloral hydrate 458
 – drinking water 454, 468
 – gas chromatography (PAT-GC) 39
 – GC/MS chromatograms 298
 – headspace technique 39ff.
 – irm-GC/MS 451
 – MAGIC 60 analysis 443
 – polymerisation additives 465
 – VOCs 437
 purity 18
 PURity soil chromatogram 326
 pyrazophos 480, 493, 505
 pyrene
 – capillary columns 119
 – *m/z* values 237
 – PAHs calibration 523
 – response gains 627
 – retention time 71
 pyributicarb 492
 pyridaben 493
 pyridafenthion 492
 pyridines 217
 pyrifeno-Z 491
 pyriminobac 492
 pyrogram 65
 pyrolysis 63ff.
 – acrylate copolymer layers 643
 – CNC compounds 154
 – FID 151
 – IRMS 198
 pyroquilon 491
 pyrrolidine derivatives 320

q

QED 489
 quadrex 121, 125
 quadrupole analyser 5, 171ff., 211
 quadrupole instruments 1, 231ff.
 quadrupole ion storage device (QUISTOR) 5
 quadrupole mass filter 5
 quadrupole mass spectrometers 84, 184
 quantitation

– acrylate copolymer layers 643
 – cannabis 589
 – clofibric acid 611
 – drinking water 454, 468
 – environmental components 625
 – mass spectroscopy 395, 401f.
 – multi-residue method 498
 – PAHs 525
 – PBBs 568
 – polychlorinated biphenyls 545
 – prostaglandins 596
 – pyrolysis 67
 – static headspace analysis 437
 – target compound analysis 477
 – toxicological-chemical analysis 607
 – trichothecene mycotoxins 621
 quantitation limit (LOQ)
 – BTEX 432
 – MAGIC 60 analysis 447
 – mass spectroscopy 397
 – PBDE analysis 561
 quasimolecular ion 213
 quaternary ammonium columns 8
 QuEChERS 489
 quinalphos 296, 492, 502
 quinidines 610
 quinines 610, 625
 quinoctamine 491
 quinoxifen 492
 quintozene 500

r

radical loss 386
 RDX mass spectra 386f.
 reactions
 – chemical ionisation 212
 – ECD 156
 reagent gases
 – charge exchange 219
 – NCI 224
 – protonation 218
 reagent ion capture 222
 rearrangements 206, 335, 343f.
 recombination 156
 recondensation 100
 reconstructed ion chromatogram (RIC) 294
 recording individual masses (SIM/MID) 233
 recoveries 432, 499
 reduction reactor 199
 redundancy filter 321, 327

- reference gas injection 204
 - reference masses calibration 264
 - reference spectra library search 332
 - reference substances 313
 - refractory methane 198
 - remoxiprid 610
 - reproducibility 63, 429, 452
 - residual gas analysers 5
 - residual monomers 465
 - residual solvents 51
 - residue analysis 8, 477
 - resins 13
 - resolution
 - capillary columns 130
 - chromatographic 140 ff., 150
 - MS 167 ff.
 - PAHs quantification 525
 - two dimensional gas chromatography 75 ff.
 - resolving power 129, 167
 - response factors 46, 223
 - Restek capillary columns 122, 126
 - restrictors 20
 - retardation term 145
 - retarding layers 643
 - retention gap 113, 116
 - retention index 309, 402
 - retention times
 - capillary columns 129
 - chromatograms 138, 149, 293
 - clofibric acid 613
 - fast gas chromatography 71
 - library search 322
 - metabolites 595
 - MID 247
 - multi-residue method 499
 - nitrogen/phosphorus pesticides 479
 - open split coupling 275
 - pesticides 491
 - prostaglandins 598
 - SIM analysis 236
 - retention volume 138, 150
 - Retro-Diels-Alder cleavage 346
 - reverse phase chromatography (RPC) 8, 9
 - RFIT soil chromatogram 324 f.
 - Rhine water 51
 - ricinine 621
 - river water analysis 532
 - rotating vane pumps 273
 - rotational perforator 505
 - RP-C¹⁸-material 611
 - rubber 409
 - run time effects 72
- s**
- salting out 35
 - sample characteristics 112, 118
 - sample injection *see* injection
 - sample inlet systems 89 ff.
 - sample loop 37
 - sample preparation 7 ff., 32, 404
 - sample quantity 54, 68
 - sample vessels 53
 - sampling 57 ff.
 - air 545
 - SPME 13
 - sarin 392
 - scan functions 232, 256
 - scan inherent calibration 247
 - scan modes 433, 445
 - scan rates 395
 - scanning device 5
 - schematics
 - 2DGC-C-IRMS instrument 540
 - BSIA/CSIA isotope analysis 195
 - cryofocusing unit 116
 - ECD 155
 - ELCD 160
 - FID 152
 - FPD 162
 - hot split/splitless injector 97
 - HRMS MID calibration 248
 - injection headspace technique 38
 - ion trap analyser 180
 - magnetic sector mass spectrometer 175
 - MCSS device 89
 - multidimensional gas chromatography 82
 - NPD 153
 - open split coupling 274
 - PDD 163
 - PFE 27
 - phenoxycarboxylic acid analysis 612
 - PID 158
 - Pt foil pyrolyser 64
 - PTV 101
 - purge and trap systems 41
 - quadrupole mass spectrometer 178
 - TOF analyser 186
 - two dimensional gas chromatography 80
 - ultra fast chromatography 76 ff.
 - schradan 501

- screening
 - confirmation analysis 548
 - environmental components 625
 - toxicological-chemical analysis 607
- sealants 545
- sealing mechanism 94
- search for identical/similar compounds (SISCOM) 320
- secbumeton 501
- sector mass spectrometers 84
- seepage water 433
- selected ion monitoring (SIM) 181, 231, 234
 - air 421
 - clofibric acid 611
 - MID 294
 - Rhine water 534
 - target compound analysis 477
- selected reaction monitoring (SRM) 253, 259, 489, 496
- selective detection techniques 1
- selective ion storage (SIS) 237, 236
- selective mass instability 6
- selectivity
 - chromatograms 144
 - ELCD 160
 - H-SRM mode 496
 - PID 160
 - stationary phases 128
 - steroid hormones 592
 - target compound analysis 477
- semivolatiles 28
- sensitivity
 - BTEX 429
 - drinking water 468
 - headspace technique 34, 53
 - mass spectroscopy 399
 - SIM analysis 233 ff.
 - target compound analysis 477
- separation
 - chromatograms 150
 - fast GC analysis 554
 - funnels 137
 - GC/MS 277
 - SPE 8
 - substance mixtures 2
- Sephadex columns 9
- septum caps bleeding 53
- septum-equipped programmable injector (SPI) 104
- sequential pyrolysis 65
- serum samples 590, 607
- sewage sludge sample analysis 559 f.
- SGE capillary columns 122, 127
- shaking frequencies 36
- signal/noise ratio 90
 - dichloroethylene 442
 - headspace GC/MS 430
 - mass chromatograms 295
 - military waste 576
 - MS/MS 251, 400
 - PAHs analysis 521
 - PTV 105
 - quadrupole analysers 173, 272
 - steroid hormones 592
 - target compound analysis 477
- significance weighting 321, 328
- silanised glass wool 104
- silica gel columns 9 f., 44
- silicon
 - capillary columns 119
 - column bleed 408
 - contaminants 414
 - isotope pattern 338, 341
- siloxane-carborane 134
- simazine
 - acquisition parameters 479
 - mass spectra 368
 - multicomponent pesticide analysis 491
 - retention times 501
- simetryn 501
- similarity match values 188
- simultaneous determination, volatile halogenated hydrocarbons 433 f.
- SISCOM (library search for identical/similar compounds) 320,
- size exclusion chromatography (SEC) 9, 526
- smoke analysis 525
- sniffing devices 89
- soft drinks 625
- soil
 - air field analysis 461
 - chromatogram 322
 - MAGIC 60 analysis 446
 - PAHs analysis 521, 582
 - static headspace analysis 437
- solid phase extraction (SPE)
 - cannabis detection 589
 - clofibric acid 611
 - drinking water 467
 - multi-residue method 498

- PAHs 521, 527
- polycyclic musks 616
- prostaglandins 596
- sample preparation 8ff.
- water 508ff.
- solid phase microextraction (SPME) 3, 12f., 26, 534
 - PBBs 568
 - plungers 13
 - VOCs 437
- solubility 16f., 118
- solution preparation (internal standard) 403
- solvent-free extraction techniques 12f.
- solvents 8ff.
 - boiling points 107
 - expansion volumes 110
 - PFE 26
 - PTV 105
 - sample injection 95ff., 103
- soman 392
- sonication 27
- sovol 541
- soxhlet
 - extraction 15, 27
 - military waste 572
 - PAHs analysis 521
 - PCDD 25
- spargers 40
- spectra (pyrolysis) 67
 - *see also* mass spectra
- spectral subtraction 298f., 302
- spectrometers 231
- spectrum display 335
- spiro molecular pumps 273
- split injection 62, 92, 96, 104f.
- splitless injection 92, 97
 - polycyclic musks 617
 - purge and trap/GC/MS systems 47
- squalene 412
- stainless steel 409, 422
- standard addition procedure 405
- standard deviations 46
- standard mean ocean water 205
- standardisation procedures 401
- static headspace chromatography (HSGC) 31
- static headspace technique 31, 49, 52
 - BTEX 429
 - volatile priority pollutants 437
- stationary phases 118, 134
- statistical errors 90
- steel plunger 13
- stepper motor 89
- steroid hormones analysis 592
- stir bar sorptive extraction (SBSE) 15
- structural isomers 313
- structure postulate 336
- structure selective detection 259
- structures, dinitrophenol herbicides 509
- styrene
 - detection limits 441
 - headspace technique 36
 - MAGIC 60 analysis 448
 - mass spectra 357
- styrol 425
- substance classes 351
- substance exchange/transport 142
- substance identification chromatograms 297
- substance transport limits 17
- substituted alkylamines/alcohols 216
- substituted phenols 472
- substituted pyridines 217
- subtraction 298f.
- sulfallate 500
- sulfides 215
- sulfotep 500
- sulfur
 - compound class 628
 - contaminants 411
 - FPD 161
 - isotope pattern 338, 341
 - proton affinities 216
- sulprofos 504
- SUMMA canister 421f.
- supelco 123, 127
- supercritical fluid extraction (SFE) 3, 8, 15ff., 582
- surface model, adsorbents 55
- surface water
 - chloral hydrate 458
 - clofibric acid 613
 - herbicides 508
 - polycyclic musks 616
- surrogate standard 611
- swelling 21
- switching 229
- sympathomimetics 603
- syringe burette sampling 59

syringe injection 5, 13
 – headspace technique 37
 – on-column 112
 syringe needle 102
 syringe pump 18
 systematic toxicological-chemical
 analysis 607

t

T2-toxin 621
 tabun 392
 tandem mass spectrometry *see* MS/MS
 target compound analysis 477
 TCDD traces 182, 224, 249
 TCMTB 503
 tebuconazole 492
 tebufenpyrad 492
 tebutan 479
 tecnacene 491, 500
 TEF 548
 Teflon (PTFE) 409
 temperature conversion elemental
 analysers (TC/EA) 194
 temperature rise times (TRT) 64f.
 temperature/time profile (TTP) 65
 temperature-programmed pyrolysis 65
 temperatures
 – BTEX determination 433
 – cryofocusing 117
 – PTV split injection 105
 – purge and trap systems 44
 – thermal extraction 68
 templated resins 13, 14
 Tenax 24, 55
 – calibration unit 61
 – BTEX determination 433
 – polychlorinated biphenyls 545
 – purge and trap systems 43f.
 TEQ 548, 554
 terbacil 480, 502
 terbufos 500
 terbumeton 479, 500
 terbuthylazine 501
 terbutryn 480, 501
 terbuthylazine 369, 479
 terpenes 309, 628
 tetrachlorobenzyltoluene 181
 tetrachlorobiphenyl 379
 tetrachlorodibenzodioxin 169, 382
 tetrachlorodibenzofuran 382
 tetrachloroethane
 – detection limits 441
 – MAGIC 60 analysis 448
 – partition coefficients 34
 – purge and trap/GC/MS systems 51
 – recoveries 432
 tetrachloroethylene
 – breakthrough volumes 57
 – irm-GC/MS 452
 – ISTD 404
 – MAGIC 60 analysis 448
 – mass spectra 354
 – purge and trap/GC/MS systems 51
 – recoveries 432
 tetrachloromethane 452
 tetrachloromethoxybiphenyl 181
 tetrachlorvinphos 503
 tetradifon 493, 504
 tetrahydrocannabinol (THC) 583ff., 589
 tetrahydrofuran 51
 tetramethrin 504
 tetrasul 503
 tetrazepam 608ff.
 theophylline 610
 thermal conductivity detector 5
 thermal decomposition 112ff.
 thermal extraction 63ff., 68
 thermal stability
 – adsorption material 54
 – injector system 112
 – mass spectra 370
 thermodesorption 3, 54ff.
 – calibration 59ff.
 – polychlorinated biphenyls 545
 – sampling tubes 58
 – soil air 461
 thermodynamic equilibrium 32
 thermolabile compounds 15, 18
 thermoscientific capillary columns 123,
 127
 thermostating 53, 59
 thick film column 421f.
 thifluzamide 492
 thiobencarb 501
 thiophene 441
 thiophosphoric acid esters 371
 Thompson mass spectrometer 4
 three-dimensional data field 293
 tilidine 610
 time/response data 85

- time-of-flight (TOF) analysers 174, 183 ff.
 - MDGC 84
 - MS 168
 - TCDD isotope pattern 184
 - velocitron 4
- TMS derivatives (opiates) 587
- tolclofos-methyl 212, 491
- tolfenpyrad 493
- toluene
 - breakthrough volumes 57
 - detection limits 441
 - environmental components 626
 - headspace GC/MS 429
 - internal standard 403
 - ionisation potentials 158, 207
 - irm-GC/MS 452
 - ISTD 404
 - m/z values 237
 - MAGIC 60 analysis 448
 - mass spectra 355, 358 f.
 - military waste 572
 - PAHs analysis 523
 - partition coefficients 34
 - PFE 27
 - purge and trap/GC/MS systems 51
 - recoveries 432
 - SFE 21
 - surface water 460
- toluol 425
- tolylfluanid 503
- total ion current chromatograms (TIC) 294, 546
- total petroleum hydrocarbons (TPH) 28
- total sample transfer 97, 102 ff.
- Toxaphen 230
- toxicological-chemical analysis 607
- toxicology library 309
- toxins quantification 621
- trace analysis 7
 - benzene hexachloride 532
 - environmental components 625
 - injector system 112
 - PBBs 568
- trace metals 9
- tralomethrin 493
- tramadol 610
- transfer times 98
- traps 44, 84
 - see also* purge and trap systems
- Trennzahl number 150
- triadimefon 312, 491, 502
- triadimenol 480, 503
- triallates 500
- triazines 368
- triazophos 492, 504
- tribufos 503
- tributyl phosphate 412
- trichloroaniline 539
- trichlorobenzenes 55, 441, 448
- trichlorobiphenyl 379
- trichloroethane
 - breakthrough volumes 57
 - detection limits 441
 - ISTD 404
 - MAGIC 60 analysis 448
 - purge and trap/GC/MS systems 51
 - recoveries 432
- trichloroethylene
 - breakthrough volumes 57
 - calibration unit 61
 - irm-GC/MS 452
 - ISTD 404
 - MAGIC 60 analysis 448
 - mass spectra 353
 - purge and trap/GC/MS systems 51
 - recoveries 432
- trichlorofluoromethane 405, 448
- trichloriodobenzene 539
- trichloropropane 448
- trichothecene mycotoxins 621
- tricresyl phosphate 415
- triethylamine 51
- trifloxystrobin 492
- trifluoromethane 432
- trifluorotrichloroethane 353
- trifluralin 500
- trihexyphenidyl 610
- trimethylbenzene 405, 448
- trimethylnaphthalene 523, 627
- trimipramine 610
- trinitroaromatics 625, 628
- trinitrotoluene (TNT) 388, 572
- trioxan 51
- triphenylphosphines 414
- triple-quadrupole technique 6, 255
- Triton-X-100 14
- tube desorption 62
- turbo molecular pumps 271
- two dimensional gas chromatography 75 ff.

u

Ugilec T 241
 ultra fast chromatography 74f.
 uniconazole 492
 unit mass resolution 178
 universal mass spectral libraries 314
 unsaturated amines/anilines 216
 urine 586, 594ff.

v

vacuum systems 270f.
 van Deemter curves 131f., 148
 vaporiser 97
 vapour pressure 16
 Varian/Chrompak 123, 127
 vegetables 498
 velocitron 4
 vernolate 499
 verrucarol 622
 vinclozolin 501
 vinyl chloride
 – breakthrough volumes 57
 – drinking water 454
 – MAGIC 60 analysis 443, 448
 – mass spectra 353, 453
 viscosity 17, 166
 Viton seal 94
 Vocarb 55
 – calibration unit 61
 – chloral hydrate 458
 – drinking water 454
 – purge and trap systems 43ff.
 volatile halogenated hydrocarbons 50
 – air analysis 421ff.
 – internal standard 405
 – MAGIC 60 analysis 443
 – mass spectroscopy 351
 – simultaneous BTEX determination 433f.
 – soil air 461
 volatile organic compounds (VOCs) 54f.
 – air analysis 421ff.
 – headspace technique 39
 – irm-GC/MS 451
 – MAGIC 60 analysis 443
 – static headspace analysis 437
 volatile phosphoric acid esters 391
 volatile priority pollutants 437
 volatile substances 69, 112
 volatiles-in-food libraries 318
 volatility 2, 7
 VPDB/VSMOW references 205
 vulcanization facilitators 409

w

warfare agents 371, 391f.
 waste water
 – clofibric acid 615
 – PBBs 571
 – polycyclic musks 616
 water 199
 – air analysis 421f.
 – chloral hydrate 458
 – cleavage 345
 – clofibric acid 611
 – eluents 8
 – herbicides 505ff.
 – MAGIC 60 analysis 450
 – NCI 226
 – on-line SPE-GC/MS 532
 – PAHs analysis 521
 – PFE 27
 – polycyclic musks 616
 – protonation 218
 – purge and trap systems 47
 – SFE 21
 – SPME 534
 – static headspace analysis 437
 – trace analysis 625
 water-soluble samples 8
 waxes 629
 weighting, library search 321, 328
 WHO 548, 554
 wide pore columns 9
 Wiley Registry of Mass Spectral Data
 316, 328

x

XAD resins 55, 61
 XMC 491
 xylene
 – breakthrough volumes 57
 – detection limits 441
 – headspace GC/MS 429
 – ionisation potentials 158
 – ISTD 404
 – mass spectra 356f.
 – partition coefficients 34
 – purge and trap/GC/MS systems 51
 – recoveries 432
 xylol 425

z

Zotepine 610